



MULTIPOINT™ TERMINAL UNITS

WALL PENDANTS M-GB

FEATURES

- Supply of medical gases, compressed air and vacuum
- Supply of electricity, low-voltage and audio-video signal
- Easy fitting of control and life support devices on shelves, rails and I.V. Poles



- Special four-sections anodized aluminum profile 250 x 116 mm, length depends from the project
- The front panels and side covers are from anodized aluminum or coated with powder paint in RAL colors
- Two railings 10 x 25 mm for accessories
- Indirect lighting, patient's zone lighting and night lighting
- Medical gas outlets, electric sockets, information signal sockets, medical gas pressure control indicators
- The panel can be equipped with different national standard outlets





SET OPTIONS:

- Shelf for the monitor
- I.V. pole for infusion , syringe pumps and drips
- Suspended shelf with drawers
- Remote control for nurse call and lighting system control
- Direct light sources for patient's examination
- Bedside screen
- Flowmeter with oxygen humidifier mounted on the rail or direct connection
- Suction systems from centralized vacuum, compressed air with fluid collection jars
- Catheter holder and other suspended devices

FIELD OF APPLICATION



- For small operating theatres
- For anesthesiology wards
- For intensive therapy wards and resuscitation rooms
- For neonatal wards
- For emergency rooms
- For treatment rooms




Wall pendants can be joined end-to-end to form a continuous module for general resuscitation rooms and intensive therapy wards



CONSTRUCTION FEATURES

- Four-section special anodized aluminum profile
- Special profile grooves are made for mounting heavy external devices
- Front covers made from aluminum profile are fixed into main profile without any fasteners
- Side covers made with pins for internal fixation
- Separate medical gas and electric sections for connection with inlets through separate sections
- Fixation of wall pendants to the walls of various constructions, medical gas connection to surface or hidden pipeline

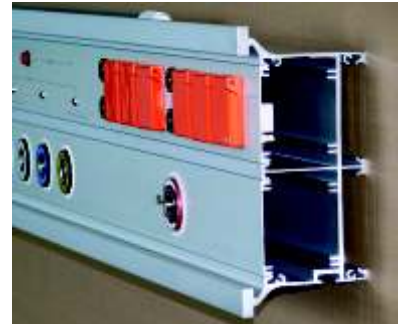

- M-GB06S Single section wall pendant for the medical gas outlets and electric sockets, with the possibility of mounting accessories and other medical



SET OPTIONS:

- Profile length 800-2300 mm
 - up to 8 medical gas outlets
 - up to 10 electrical outlets
 - up to 6 pcs equipotential outlets
 - up to 4 connectors of low voltage signal for connection of control equipment and data transferring to the information center
- Can be equipped with:
 - rails for accessories
 - suspended shelves with drawers
 - shelves for monitor and other accessories



Colors of front panels and side covers:



Color of profile: Color of shelves:



VERTICAL WALL PENDANTS



FEATURES

- For two beds service in resuscitation, intensive care wards with maximum space saving
- Supply of medical gas, electricity and data signal for each patient
- I.V. poles for infusion, syringe pumps and drips
- Upper and lower indirect lighting around the patient
- The panel can be equipped with different national standard outlets.

Colors of side panels:



● Colors of profile:



● Colors of shelves:



FLUSH-MOUNTED WALL PENDANTS



FEATURES

- Flush mounted installation
- For wards of improved design and comfort
- Supply of medical gas, electricity and data signal
- General indirect lighting, direct lighting for reading and night lighting with LED lamps
- Wireless remote for nurse call and lighting control
- Prompt setup of accessories
- Length of the profile depends from the number of gas outlets, electric sockets and sockets of data signal



● Colors of front panels:



Color of profile:



Color of shelves:



2 YEARS GUARANTEE for all wall pendants.

Wall pendants are designed and manufactured in accordance with directive 93/42/EEC requirements.

Email: mtlbi@medtech.lt